


User Manual

MS-32J4

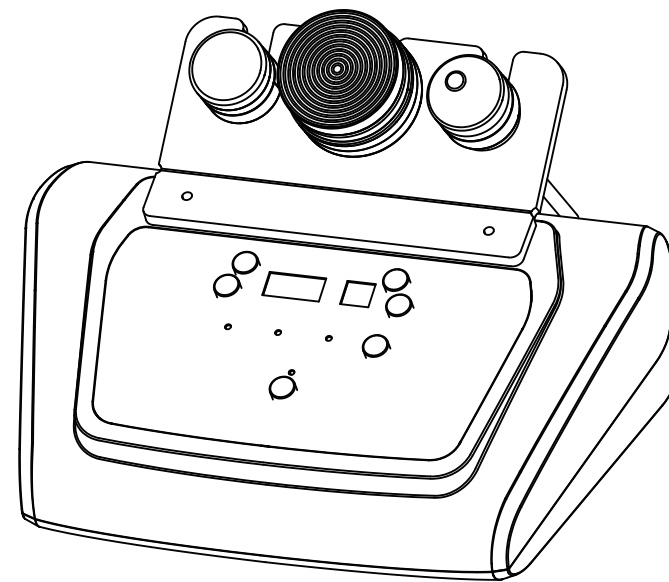

Content

1.Safety Precautions	1
2.Contraindications	1
3.Overview	2
4.General Description	3
5.Technical Specification	4
6.Use Precautions	4
7.Maintenance & Cleaning	5
8.FAQ	6
9.Ask & Answer	8
10.Operation Guide	8
11.Warranty & Support	11

1.Safety precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 7.Ensure that the instrument is applied on moist skin and avoid the operation on dry skin.
- 8.Power off the instrument when not in use to ensure its safety.
- 9.If you are currently taking anti-adiposity drugs, please discontinue their use for 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 10.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

2.Contraindications


Before using this instrument, it is essential to consult your doctor or seek advice from professionals. Details below:

- 1.Pregnancy or lactation.
- 2.Hypertension, heart disease, or with a pacemaker installed.
- 3.Unhealed wounds or during the recovery phase of surgery.
- 4.Epilepsy, severe diabetes, and hyperthyroidism.
- 5.Malignant tumors, hemophilia, and hemorrhage.
- 6.Skin disease or infectious disease.
- 7.Cesarean section within the last six months, a miscarriage within the last three months, or a natural birth within the last two months.
- 8.Had plastic surgery within 3 months or with implants or metal objects in the treatment area.





3.Overview

With the most popular and safest technology for anti-aging, firming skin, and body shaping, this 3-in-1 beauty instrument ensures you achieve your desired body figure and radiant complexion. It provides the perfect solution for people seeking improvements in skin appearance and body figures without injections, drugs, or side effects. This external device delivers discernible improvements with each passing day.

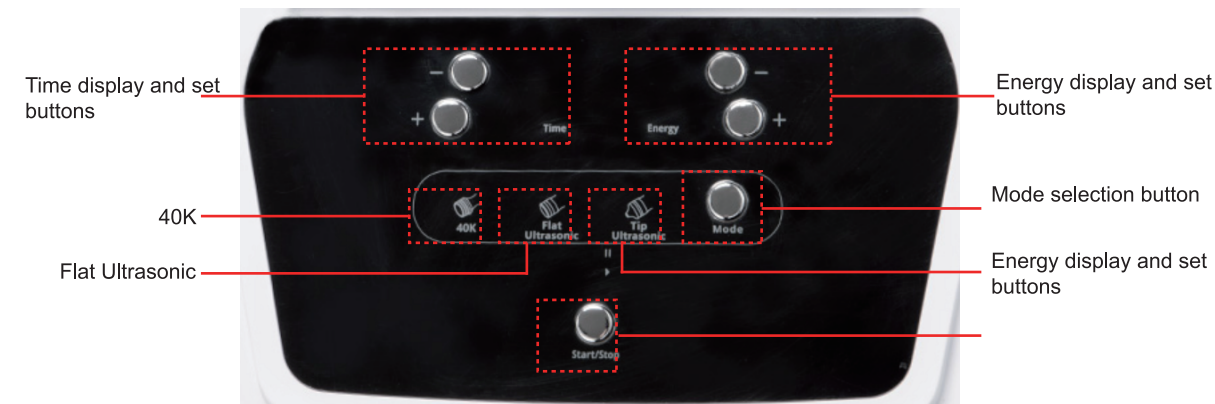4.General Description

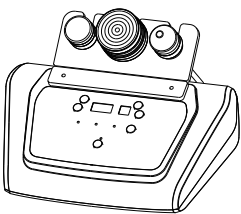

5. Technical Specification

	Power supply Input: AC110V-240v 50Hz-60Hz		Net Weight: 2.6KG		
	Power: 12W		Gross Weight: 3.2KG		
	Dimensions: 13.1×8.7×3.7 inches				
	Input Power: 12W				
	Power: 12W			Power: 12W	
	Frequency: 1MHz			Frequency: 1MHz	
				Power: 7W	
				Frequency: 40KHz	

6. Use Precautions

1. Strictly follow the instructions in the manual to use the instrument and train operator.
2. Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.

3. Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
4. Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
5. Discontinue the use of anti-adiposity drugs 1-2 months before the treatment. If immediate treatment is desired, consider extending the regimen.
6. Do not use the device on an empty stomach, and take the treatment 1 hour later after the meal.
7. Ensure full skin contact with the treatment head throughout the session to avoid discomfort.

7. Maintenance & Cleaning

1. Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
2. Clean and sanitize the applicator with a damp cloth or saline, and then dry it thoroughly with a dry towel to ensure it remains dry.
3. After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
4. Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

8. FAQ

1. Liposuction or this device?

A: Liposuction surgery utilizes negative pressure to remove adipose tissue accumulations from specific areas, resulting in a body-shaping effect. It provides you with immediate results, but it also carries risks, requires anesthesia, and involves downtime. While our device offers a risk-free solution with no side effects. It sculpts your body, stimulates collagen regeneration for lifted skin, and creates the ideal figure for an overall improved appearance.

2. Can I apply ultrasound on the entire body?

A: NO. Ultrasound disrupts adipose tissue, which is then absorbed and metabolized by the human body. The power of ultrasound is quite strong, and the heart is particularly sensitive to sound waves due to being a hollow organ. Sound waves can cause vibrations in the heart, and there can be variations in how sound waves are conducted through muscle tissue and blood within the heart, resulting in back-and-forth reflections. As a consequence, the valve may become dislodged from the heart muscle. Certainly, ultrasound treatment should be strictly avoided around the eye area to prevent the risk of retinal detachment, as well as areas such as the waist, back, chest, and the region around the heart.

3. Are there any side effects of ultrasonic treatment?

A: The ultrasound treatment is a non-surgical, non-invasive procedure that requires no incisions or anesthesia. It combines various technologies, including ultrasound (mechanical wave), cavitation, and ultrasound focusing. This procedure delivers multiple benefits such as breaking down adiposity, dissolving adiposity, tightening the skin, and sculpting the body. It selectively targets low-density adipose tissue while safeguarding high-density tissues like blood vessels and nerves, ensuring no adverse effects on the body. You may notice a slight noise in your ears during the treatment, which is

a normal occurrence and requires no cause for concern.

4. Why do I get tinnitus during the ultrasonic treatment?

A: When the ultrasound with concentrated energy vibrates strongly and has a frequency higher than 20KHz, it can penetrate the skin to reach the adipose layer that is 20mm deep. The energy causes high-speed friction and heat within the adipose tissue, which breaks them down and emulsifies them. However, this process can create noise that may be audible in your ears.

9. Ask & Answer


Ask	Cause	Answer
The 40K handle is not functioning.	You might mistakenly assume that the handle is not functioning properly because of its low-noise design.	To check if the handle is functioning, please refer to the 40K handle test.
	You may have forgotten to apply the gel or have used a product that is not designed for use with the instrument.	Please use the specific gel designed for the 40K treatment.
	Insufficient pressure.	Apply more pressure during 40K handle application.

10. Operation Guide

Face:



Step 1: (10 minutes)

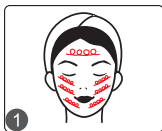

Eye:




Step 1: (10 minutes)



Arm:




Step 1: (15 minutes)

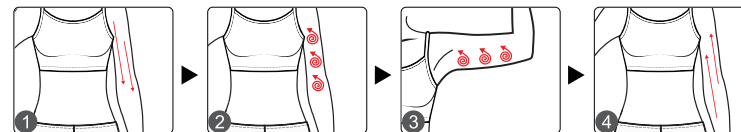

Belly:




Step 1: (15 minutes)

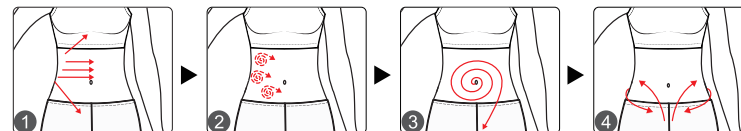

Leg:



Step 1: (15 minutes)

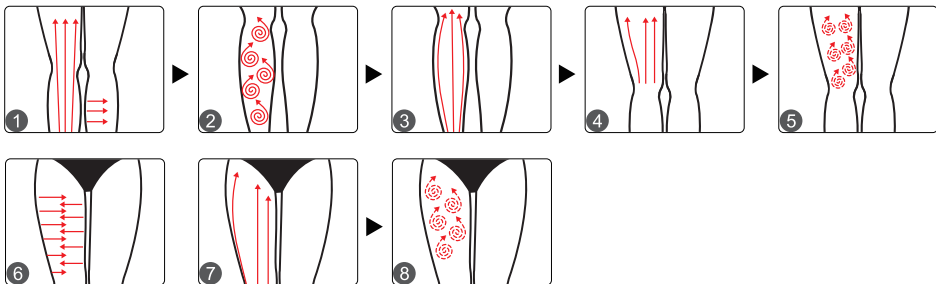

11. Warranty & Support

We provide a one-year warranty and a one-month return and replacement policy from the date of purchase.

Product Name: Versatile Beauty Enhancer
 Product Model: MS-32J4
 Production Date: _____
 Quality Inspection: _____

PASS